




Installation/Service Manual
Circulating Water Baths
Models 260 (2864/2865), 265 (2866/2867), 270 (2868/2869)

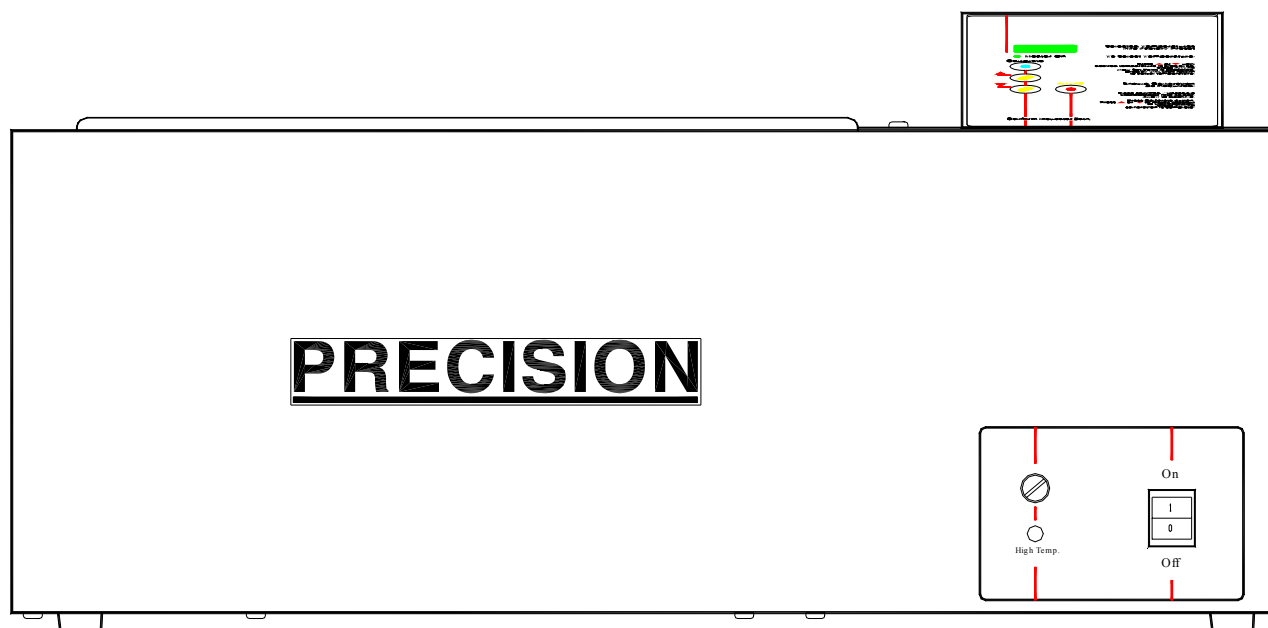

Thermo Electron Corporation
Millcreek Road, Box 649
Marietta, Ohio 45750
USA
Phone: 740-373-4763
Toll Free: 800-848-3080
FAX: 740-373-4763

Manual P/N 3177890
Rev. E Dated 18 April, 2005

NOTICE

THE MATERIAL IN THIS MANUAL IS FOR INFORMATION PURPOSES ONLY. THE CONTENTS AND THE PRODUCT IT DESCRIBES ARE SUBJECT TO CHANGE WITHOUT NOTICE. THERMO ELECTRON CORPORATION MAKES NO REPRESENTATIONS OR WARRANTIES WITH RESPECT TO THIS MANUAL. IN NO EVENT SHALL THERMO BE LIABLE FOR ANY DAMAGES, DIRECT OR INCIDENTAL, ARISING OUT OF OR RELATED TO THE USE OF THIS MANUAL.

For repair information or replacement parts assistance from the manufacturer, call Technical Services using our toll free telephone number.

888-213-1790
(FAX) 740-373-4189

REVISION STATUS

INDEX	DATE	AMENDED PAGES	NOTES
A	6/99		Initial release
B	5/00	16	Update pump assembly kit numbers (new pumps)
C	NOV01	6	Add caution "acidic & caustic substance"
D	APR02	5,8,27	Add inlet screen info, correct wire color
E	APR05	36100125 (34002503)	new manual #, manufacture location

Contents

INTRODUCTION	1
UNPACKING AND DAMAGE	1
GENERAL INFORMATION.....	2
PERFORMANCE	2
INSTALLATION	2
WATER LEVEL REGULATOR	3
EXPLANATION OF CONTROLS	6
OPERATION	7
MAINTENANCE.....	8
TROUBLESHOOTING	8
PARTS REPLACEMENT	13
TEMPERATURE CALIBRATION	16
REPLACEMENT PARTS LIST	17
ASSEMBLY & SCHEMATIC DWGS	19
WARRANTY	28

INTRODUCTION

Your satisfaction and safety are important to **Thermo** and a complete understanding of this unit is necessary to attain these objectives.

As the user of this apparatus, you have the responsibility to understand the proper function and operational characteristics of your bath. This instruction manual should be thoroughly read and all operators given adequate training before attempting to place this unit in service. Awareness of the stated cautions and warnings, and compliance with recommended operating parameters — together with maintenance requirements — are important for safe and satisfactory operation. The unit should be used for its intended application; alterations or modifications will **VOID THE WARRANTY.**

WARNING

*AS A ROUTINE LABORATORY PRECAUTION,
ALWAYS WEAR SAFETY GLASSES WHEN
WORKING WITH THIS APPARATUS.*

This product is not intended, nor can it be used, as a sterile or patient connected device. In addition, this apparatus is not designed for use in Class I, II or III locations as defined by the National Electrical Code.

UNPACKING AND DAMAGE

Save all packing material until unit is put into service. This merchandise was carefully packed and thoroughly inspected before leaving our factory.

Responsibility for safe delivery was assumed by the carrier upon acceptance of the shipment; therefore, claims for loss or damage sustained in transit must be made upon the carrier by the recipient as follows:

1. **Visible Loss or Damage:** Note any external evidence of loss or damage on the freight bill, or express receipt, and have it signed by the carrier's agent. Failure to adequately describe such external evidence of loss or damage may result in the carrier's refusing to honor your damage claim. The form required to file such claim will be supplied by the carrier.
2. **Concealed Loss or Damage:** Concealed loss or damage means loss or damage which does not become apparent until the merchandise has been unpacked and inspected. Should either occur, make a written request for inspection by carrier's agent within fifteen (15) days of the delivery date; then file a claim with the carrier since the damage is the carrier's responsibility.

If you follow the above instructions carefully, we will guarantee our full support of your claim to be compensated for loss or concealed damage.

DO NOT — FOR ANY REASON — RETURN THIS UNIT WITHOUT FIRST OBTAINING AUTHORIZATION. In any correspondence to **Thermo**, please supply the nameplate data, including catalog number and serial number.

GENERAL INFORMATION

Precision Circulating Baths are widely used in research and quality control. Their superb temperature uniformity and stability makes them especially desirable for legal or reference tests.

The microprocessor control panel houses all functions necessary to operate the bath. The 5 push-button switches and single digital display window allow the operator to adjust bath temperature and temperature calibration via a single set of controls.

The proportional integral temperature control coupled with agitation supplied by a centrifugal pump located in the control housing allows precise temperature control. Use of the gable cover provided is required to maintain optimal temperature sensitivity.

A high limit thermostat is provided and can be set to prevent heater runaway in the event of temperature control failure.

The interior of the bath is constructed of stainless steel and is designed for operation with distilled water or water solutions, such as water /ethylene glycol with corrosion inhibitors added. The body is made from galvanized steel and is painted for added protection. A drain is located at the far left hand end of the bath. A stainless steel gable cover is also provided with the bath.

The 230-volt units are identical in appearance to the 115-volt units.

PERFORMANCE DATA

The following table identifies the specifications for the Circulating Baths.

4.02 Performance Data			
	With Cover		
	Model #260	Model #265	Model #270
Uniformity 37°C	± .05	± .05°C	± .05
Sensitivity 37°C	± .05	± .05°C	± .05

INSTALLATION

WARNING

INSTALLATION SHOULD BE COMPLETED BY QUALIFIED PERSONNEL ONLY.

Location - The most uniform operating conditions will be obtained by placing the bath on a level surface in an area remote from drafts, ventilating outlets, radiators, and other rapidly changing ambient conditions.

Electrical Connections -

IMPORTANT

For personal safety, this apparatus must be properly grounded.

Table 1 Listing of Models Included In This Manual					
Category Number	Model	Volts	Hertz	Watts	Amps
3166680	#260	120	50/60	1000	8.3
3166682		230	50/60	1000	4.3
3166684	#265	120	50/60	1550	12.9
3166685		230	50/60	1550	6.7
3166681	#270	120	50/60	1550	12.9
3166683		230	50/60	1550	6.7

1. The power cord provided on this unit is equipped with a three-prong (grounding) plug which mates with standard three-prong grounding wall receptacle to minimize the possibility of electric shock hazard from this apparatus. If in doubt the user should have the wall receptacle and circuit checked by a qualified electrician to make sure the receptacle can provide adequate current and is properly grounded.
2. Where a standard two-prong wall receptacle is encountered, it is the personal responsibility and obligation of the user to have it replaced with a properly grounded three-prong wall receptacle. *Do not, under any circumstances, cut or remove the third (ground) prong from the power cord. Do not use a two-prong adapter plug.*

Determine the total amount of current being used by other apparatus connected to the circuit that will be used for this apparatus. It is critical that the added current demand (see nameplate) of this and other equipment used on the same circuit does not exceed the rating of the fuse or circuit breaker.

CAUTION

BE SURE THAT THE POWER SUPPLY IS OF THE SAME VOLTAGE AS SPECIFIED ON THE NAMEPLATE.

TO AVOID BUILD-UP OF MINERAL DEPOSITS AND TO PREVENT CORROSION, USE ONLY DISTILLED WATER (50K OHMS TO 1-M) IN THE BATH. DO NOT USE DEIONIZED WATER, OVERTIME IT BECOMES CORROSIVE AND WILL ATTACK STAINLESS STEEL. IF DEIONIZED WATER IS THE ONLY TYPE AVAILABLE THEN USE A MIXTURE OF WATER WITH A RATIO OF 50% DEIONIZED WATER AND 50% TAP WATER.

ELECTROLYSIS CAN DAMAGE STAINLESS STEEL. THIS CAN OCCUR IF AN OBJECT IS ALLOWED TO REST DIRECTLY ON THE SURFACE, TRAPPING MOISTURE THAT BECOMES OXYGEN STARVED BUT IS SURROUNDED BY WATER CONTAINING OXYGEN. THE RESULTING ELECTROLYTIC ACTION WILL PIT OR ERODE THE STAINLESS STEEL.

WATER LEVEL REGULATOR **(OPTIONAL)**

INSTALLATION INSTRUCTIONS - Use the following steps along with Fig. 1 & 2 shown on the following page.

1. Turn off unit, disconnect power & drain water from the bath.
2. Unscrew drain plug from side of bath body.
3. Replace with elbow connector (shipped loose).

NOTE


TO INSURE GOOD SEAL, WRAP THREADS WITH TEFLON TAPE OR EQUIVALENT SEALER.

4. Tighten elbow connector so open end (tube fitting) faces straight down.
5. Place Water Level Regulator Bkt. Ass'y on the lip of the bath pan (drain side of the bath towards the rear of the unit). The small flange with gasketing near the top of the bracket should rest against the inside of the bath pan. The other 2 flanges with gasketing should rest against the outside of the bath body.
6. Insert the plastic tubing (coming from the bottom of the water level regulator) into the elbow connector (at the drain) & tighten the elbow connector fitting firmly to insure a good seal.
7. Uncoil Tygon tubing & press onto 1/4" Hose Barb Fitting (located at the top of water level regulator).
8. Hook-up other end of the tygon tubing to water supply & turn on water. (Water pressure should fall within the following range: MIN is approx. 15 PSI, MAX is approx 75 PSI)


WATER LEVEL REGULATOR ASSEMBLY

NOTE
WATER LINE MARKING ON THE WATER LEVEL REGULATOR IS AN INDICATOR OF THE APPROX. WATER LEVEL THAT WILL BE MAINTAINED IN BATH. THE ASS'Y HAS BEEN SET TO MAINTAIN APPROX. 1-3/4" OF WATER. TO RAISE THE SET WATER LEVEL, LOOSEN THUMBSCREW & RAISE WATER LEVEL REGULATOR TO DESIRED LEVEL (COINCIDE WITH WATER LINE MARK) & TIGHTEN THUMBSCREW TO HOLD WATER LEVEL REGULATOR IN PLACE.

NOTE:
DO NOT OVER-TIGHTEN THUMBSCREW

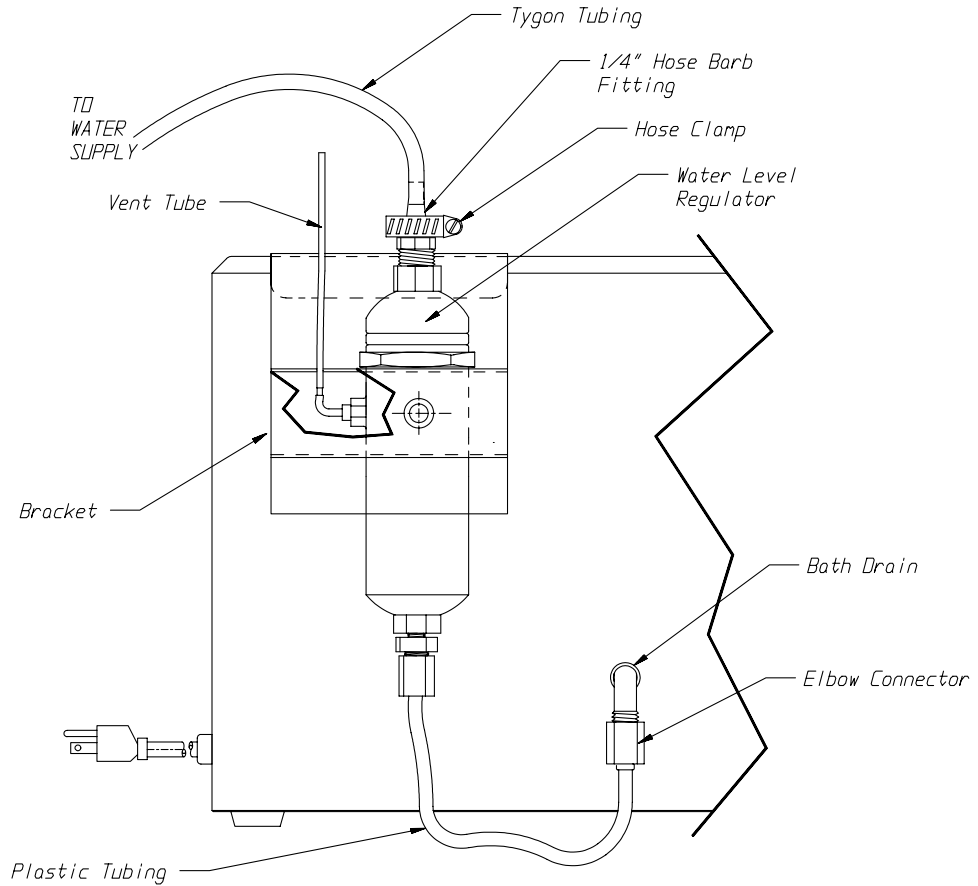

**FIGURE 1
SIDE VIEW OF BATH**

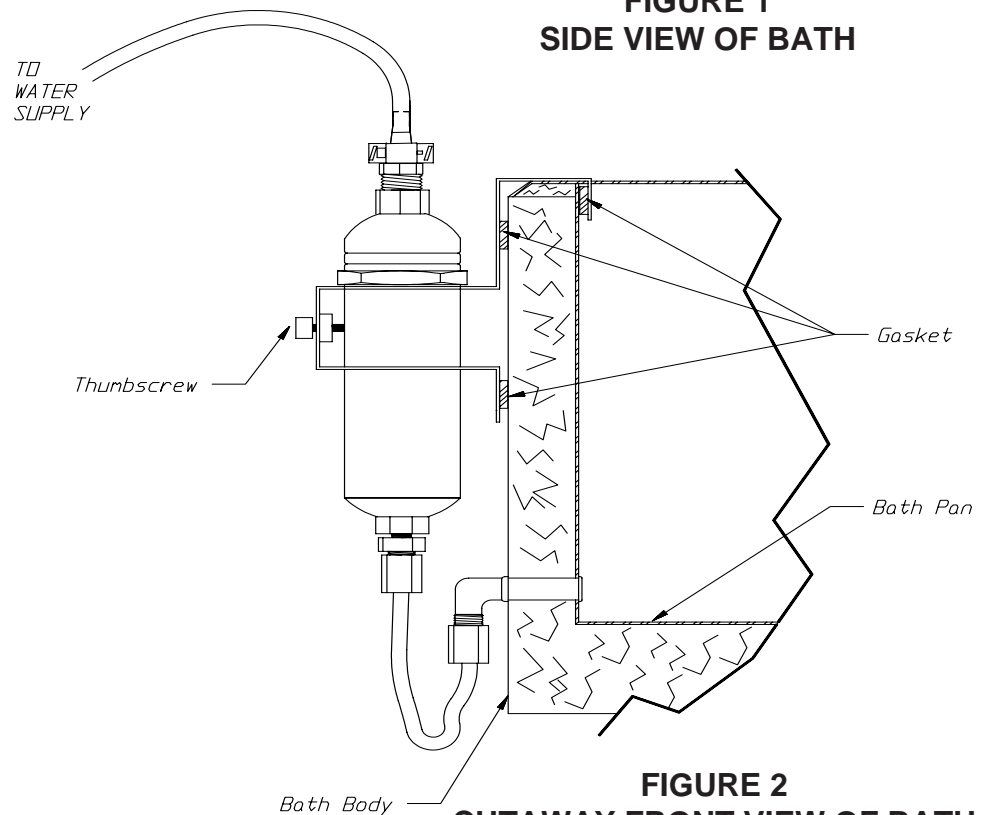


**FIGURE 2
CUTAWAY FRONT VIEW OF BATH**

INLET SCREEN ASSEMBLY **(OPTIONAL)**

An optional inlet screen assembly is packaged with each unit. If used, this screen will filter foreign materials (i.e. pieces of paper) from being sucked through the pump. To install the screen assembly, refer to the diagram below.

NOTE:

*IF USED, THE SCREEN ASSEMBLY MUST
BE INSPECTED DAILY. REMOVE ANY
DEBRIS FOUND ON THE SCREEN.*

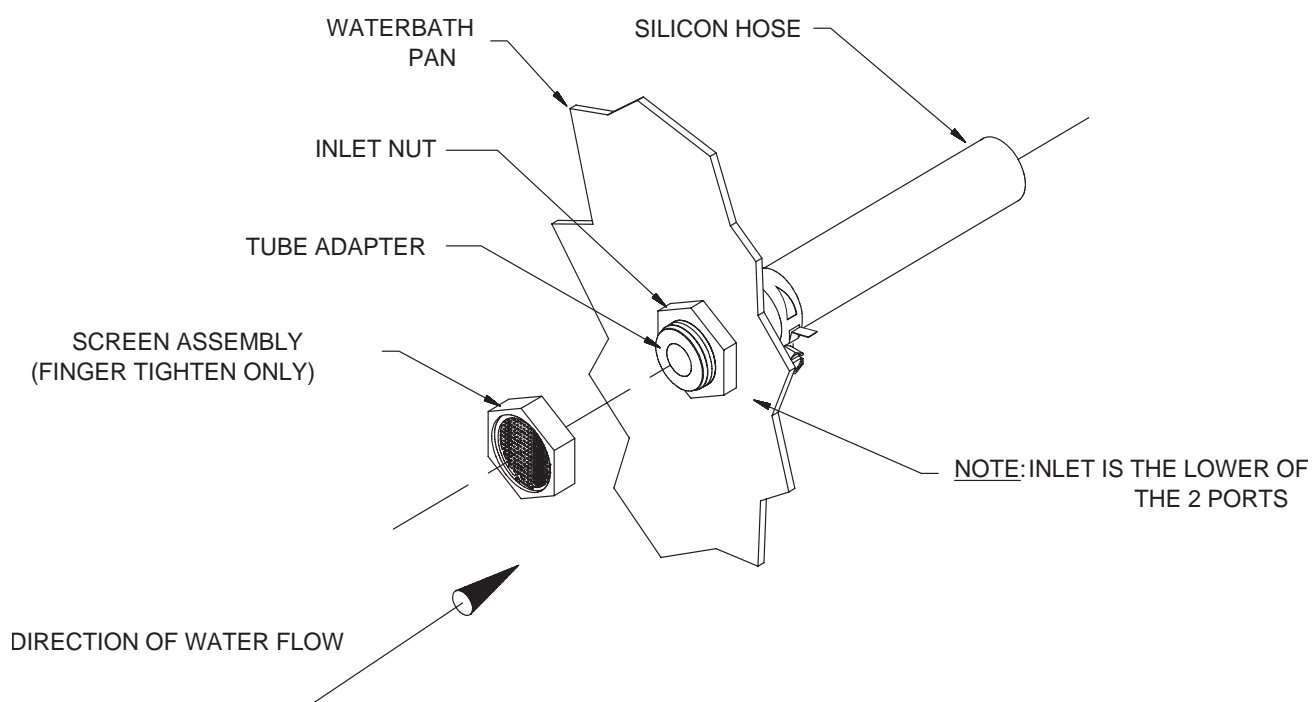

EXPLANATION OF CONTROLS

Power Switch - The power switch is located on the lower right hand side of the bath, it provides power for the entire unit.

Hi-Limit Thermostat - is located immediately to the left of the power switch on the lower right hand corner of the bath. The high limit is an adjustable control used to prevent thermal runaway in the event that the primary control fails. After you have stabilized the bath temperature, use a standard screwdriver to adjust the potentiometer. When delivered, the Hi-Limit is turned fully clockwise (highest setting). In this position the Hi-limit light will be **OFF**. Turn the potentiometer counter-clockwise until the Hi-Limit light is **ON**, and then adjust clockwise 1/8 of a turn. Verify that the Hi-Limit light is **OFF**. Now if your bath temperature rises above the desired setting the Hi-Limit protection will be activated.

Control Panel

The Control panel is located on top of the bath and provides primary temperature control. The Control Panel contains the following features (See Figure A)

1. LED Display

4-digit display used to show both actual and set point values for the temperature. Also shows configuration parameters when in configuration mode.

2. Heater on Lamp - the Heater on Lamp is lighted when power is applied to the heating element.

3. Indicator Lamps - These lamps indicate which parameter (Actual Temp or Setpoint Temp) is displayed in LED window.

4. Set Point Key - The Set Point Key is used to toggle the display between actual bath temperature and the setpoint value.

5. Temperature Selection Keys


The Temperature Selection Keys are used to select the desired temperature setpoint value and are only active after the Set Point Key has been pressed.

6. Enter Key

The Enter Key is used to store a new set point value.

7. Offset Key

The Offset Key is used to adjust the actual temperature value to match a calibrated reference thermometer.

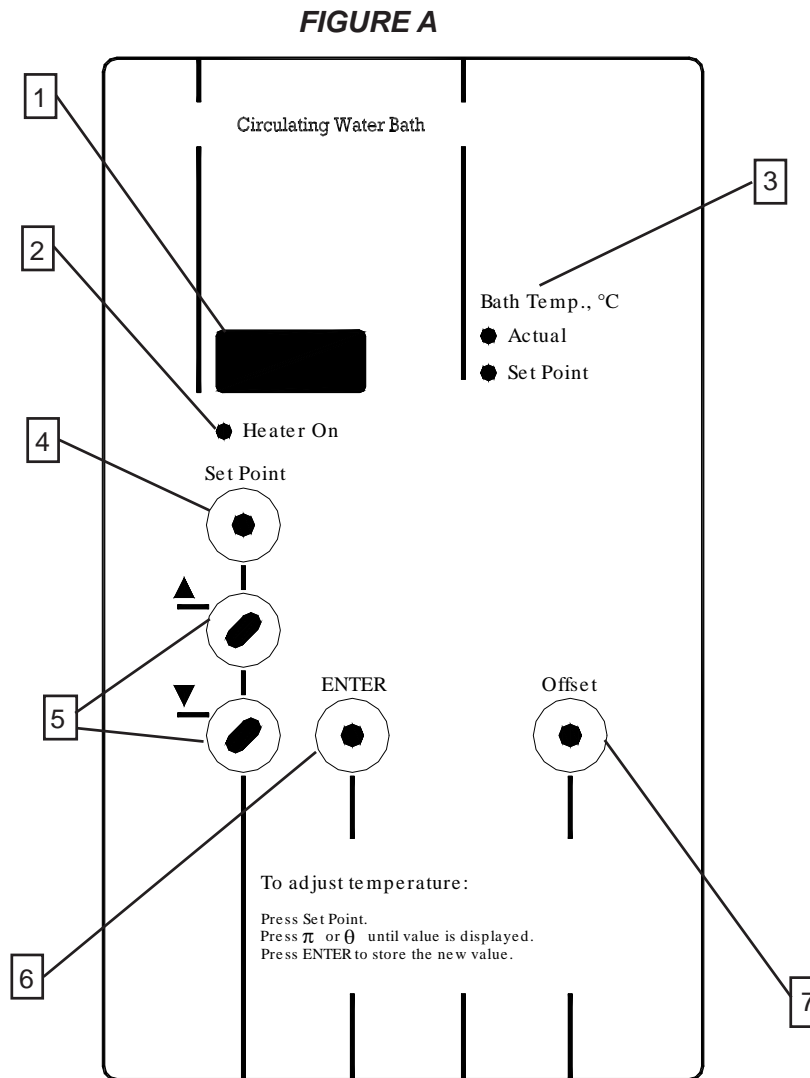

OPERATION

CAUTION:

EXERCISE CARE WHEN USING ACIDIC OR CAUSTIC SOLUTIONS AS THEY WILL ATTACK THE GALVANIZED STEEL BATH BODY IF SPILLED INTO THE BATH. IF SPILLS DO OCCUR, THE BATH LIQUID SHOULD BE IMMEDIATELY DRAINED AND THE UNIT THOROUGHLY FLUSHED. SPILLS AND CONDENSATION SHOULD BE CLEANED/ REMOVED FROM ALL METAL SURFACES AFTER EACH USE.

When filling the bath with water, allowance must be made for the displacement of water upon immersion of samples.

CAUTION

WATER LEVEL REQUIREMENTS

MINIMUM LEVEL - WATER LEVEL SHOULD BE 1 1/2 INCHES ABOVE THE PUMP OUTLET. THIS WILL KEEP THE PUMP FROM RUNNING DRY AND FAILING TO OPERATE.

MAXIMUM LEVEL - WATER LEVEL SHOULD NOT BE HIGHER THAN 1 INCH FROM THE TOP OF THE BATH PAN.

To conserve energy, reduce evaporation and increase temperature control accuracy, use the gable cover provided. Do not use aluminum foil as a cover, as it may cause corrosion due to an electrochemical reaction.

1. Add distilled water to bath. Use water levels as instructed above.
2. Power up: Depress the power switch located at the lower right hand corner of the bath. The unit is fully operational upon application of power and the display will read actual bath temperature
3. **To set temperature control point:**
 - a. Press SETPOINT. The display will now show the current temperature setpoint.
 - b. Press UP or DOWN until the desired temperature setpoint is displayed.
 - c. Press ENTER to store the new value.

- d. Display will return to actual temperature.
- e. If ENTER is not pressed, new value will not be set and display will return to actual temperature.

4. To adjust temperature offset:

- a. Temperature offset should be used to match a stable bath's actual temperature display to a calibrated thermometer.
- b. Press ENTER until actual temperature is displayed.
- c. While pressing OFFSET, press UP or DOWN until display shows the **DIFFERENCE** between displayed temperature and reference thermometer. Continue holding offset and press ENTER to set new value.
- d. Display should now match your calibrated thermometer.
- e. If ENTER is not pressed, new value will not be set and display will return to previous actual temperature.
- f. Example: Setpoint is 37° and actual bath temperature display is stable at 37°. Calibrated reference thermometer reads 37.2°. Press OFFSET and UP until display reads .2. Press ENTER while still holding OFFSET key. Release keys. Actual bath temperature will read 37.2°. Readjust setpoint to 37° by repeating step 3.

NOTE

AS WATER BATHS ARE OPERATED, ESPECIALLY AT ELEVATED TEMPERATURES, CONDENSATION WILL FORM ON THE UNDERSIDE OF THE GABLE COVER. THE AMOUNT OF CONDENSATION CAN BE CONSIDERABLE. CARE SHOULD BE EXERCISED WHEN REMOVING THE COVER SO THAT WATER WILL FALL BACK INTO THE BATH CHAMBER, AND NOT ON THE BATH SURFACE. THE COVER SHOULD BE PLACED UPSIDE DOWN WHEN REMOVED TO MAINTAIN A DRY WORK AREA.

MAINTENANCE

Cleaning: Stainless steel will resist corrosion; however, it is not impervious to it. Proper maintenance of the stainless steel bath chamber will help assure many years of service.

CAUTION

AVOID SPILLING HARSH CHEMICALS ONTO THE BATH, AS CORROSION OF THE STAINLESS STEEL MAY RESULT.

It should be cleaned regularly with mild soapy water and rinsed with distilled water. Should algae or other undesirable microorganisms form on the top of the bath media, add a little formaldehyde or quaternary ammonium germicide, available from Precision Scientific (P/N 51200912).

IMPORTANT

IF IT IS NECESSARY TO USE THE FOLLOWING CHEMICALS, LIMIT THE TIME TO A MAXIMUM OF FOUR HOURS. CLEAN SURFACES IMMEDIATELY AFTER USE.

-Aluminum Chloride	-Barium Chloride
-Bichloride of Mercury	-Calcium Chloride
-Carbolic Acid	-Chlorinated Lime
-Citric Acid (boiling)	-Dakin's Solution
-Ferrous Chloride	-Mercury Salts
-LysolMercuric Chloride	-Phenol
-Potassium Permanganate	-Tartaric Acid
-Potassium Thiocyanate	-Stannous Chloride
-Sodium Hypochlorite	

CAUTION

NEVER USE THE FOLLOWING CHEMICALS:

-Aqua Regia	-Ferric Chloride
-Iodine	-Sodium Azide
-Sulfuric Acid	

Removing Discoloration: Should the stainless steel ever become discolored by iron rust, use the following procedure to remove all traces of the rust and restore the stainless steel.

WARNING

OBSERVE THE FOLLOWING SAFETY PRECAUTIONS! USE HEAVY GLOVES OR OTHER ADEQUATE HAND PROTECTION. WEAR GOGGLES OR OTHER ADEQUATE EYE PROTECTION. ONLY WORK IN AREAS WITH ADEQUATE VENTILATION.

Prepare a solution of 20% nitric and 1.5% hydrochloric acid (if preferred, a 2% to 5% solution of warm oxalic acid may be used). Swab solution over surface, allowing it to remain until all rust is loosened. This will usually take 1 to 2 minutes.

As soon as rust is loosened, immediately flush with clean water until all acid is removed.

Screen assembly option: If this option is installed, inspect daily and remove any debris from the screen.

TROUBLESHOOTING

WARNING

*SERVICE SHOULD BE PERFORMED BY A QUALIFIED TECHNICIAN. BEFORE REPLACING ANY ELECTRICAL OR MECHANICAL COMPONENTS, UNPLUG THE LINE CORD. IF ELECTRICAL POWER IS REQUIRED FOR SERVICE, **USE EXTREME CARE.***

Refer to Problems and Solutions for troubleshooting information on the baths. This guide provides the basic information required to repair the bath.

The following is a list of the tools and instruments required to perform the procedures outlined in the Troubleshooting Procedures table.

Tools Needed:

- Phillips or Flat Blade Screwdriver
- 7/16" Socket or adjustable wrench
- Ohmmeter
- DC Voltmeter
- AC Voltmeter
- 5VDC Power Supply

When changing circuit boards, please use CAUTION when re-installing the flat cable connectors that come from the display board. Make sure the pins are NOT bent or broken. Boards that are otherwise fine will not operate correctly if these pins have been abused because they will not receive the necessary signals from other boards.

Before attempting any troubleshooting for a particular problem, it is good practice to verify power supply voltages of the POWER SUPPLY board #34372501. Read these steps before performing them to get acquainted with the procedure.

1. Be sure that the main power source is OFF or disconnected.
2. Connect J101 and J105 to the POWER SUPPLY BOARD. DO NOT connect J102 that comes from the CPU and the DISPLAY at this time.
3. Connect the negative lead of the voltmeter to TP2 - COM. Select a VDC of 20. Connect the positive lead to TP1 - UNREG.
4. Turn ON the main power.
5. The voltmeter should be reading approximately 14 volts DC ± 2 .
6. Connect the positive lead to TP3 $\pm 5V$. The voltmeter should be reading +5 volts DC ± 0.5 .
7. Turn OFF main power.

8. Connect the display board ribbon connector to the power supply board at J102 and then repeat steps 3 thru 6. Do not connect the connector to the CPU board at this time. If the above voltages are not measured again, the display board has a bad component which is loading down the voltages. If the voltages are still present, continue on to step 9.
9. Connect the display board ribbon connector to the CPU board J310 and repeat steps 3 thru 6. If the above voltages are not measured again, the CPU board has a bad component which is loading down the voltages. This board will have to be replaced.

GENERAL INFORMATION: The POWER SUPPLY PCB ASSEMBLY #34372501 has the triac (solid state AC voltage switch) on it which supplies the power to the heater. This triac is "told" to operate, when need be, by the microprocessor. There is another device known as an opto-coupler which serves as the high/low voltage isolator between the triac and the microprocessor.

When the bath starts experiencing temperature problem and/or variations, due to known controller malfunctioning, these two components become prime suspects. The reason being they are under higher operating stress than other components.

If the bath starts experiencing temperature problems such as no heat, or "creeping" upward heat, then follow the troubleshooting instructions below. These instructions require the use of a voltmeter being able to measure DC and AC voltages (preferably digital) and an ohmmeter.

WARNING

THE FOLLOWING TROUBLESHOOTING INSTRUCTIONS REQUIRE THAT POWER BE ON. ONLY QUALIFIED SERVICE PERSONNEL SHOULD PERFORM THESE PROCEDURES.

Open the control cover and familiarize yourself with the POWER SUPPLY assembly #3166924. Locate the triac (Q2), the opto-isolater (U1), the resistor (R4), and the test point #2 (TP2-COM). The first measurements will be DC voltage measurements. A DC scale of at least 10 volts should be selected.

PROBLEMS & SOLUTIONS

******* Problem A : No Heat *******

1. Verify that the setpoint temperature is greater than the actual water temperature.
2. Verify that the HIGH TEMP light is "OFF". If the light is "ON", turn the HIGH LIMIT CONTROL fully clockwise.
3. Check temperature probe voltage reference. These measurements are on the CPU board #3167119.
 - a. Place the negative lead of a DC voltmeter on TP1-COM test point. Measure the voltage at U7-pin 2. It should be approximately 14 volts DC. Now measure the voltage at U7-Pin 6. It should be approximately 5 volts DC. If the 14 volts is present but not the 5 volts, the CPU board must be replaced.
4. Check temperature probe.
 - a. Measure the bath water temperature and make note of it. Make sure it will not change drastically during this check.
 - b. With a DC voltmeter, measure the supply of the temperature sensor at J302 pins 1 & 3. It should be approximately 5 volts. If it is, then continue on to step C. If not, see Step 3a.
 - c. Measure the output of the temperature sensor at J302 pins 2 & 3. Pin 3 being the reference or ground. The output relationship of voltage to degrees centigrade is listed in the table below. Compare the measured voltage to the temperatures listed to determine if the probe is working.

Temp Degrees C	Volt	Temp Degrees C	Volt
10	0.500V	60	1.400V
20	0.680V	70	1.580V
30	0.860V	80	1.760V
40	1.040V	90	1.940V
50	1.220V	99	2.218V

5. The bath is not heating when it should be requesting heat. Voltage measurements are made on the POWER SUPPLY BOARD #3166924.
 - a. Select a setpoint temperature at least 10° above what the actual temperature is. The bath should be requesting heat as indicated by the heater "ON" indicator lamp on the front panel. It should be on continuously, not flashing.
 - b. Measure the voltage at U1-PIN 2 with respect to TP2-COM. It should be no greater than 4 VDC.
 - c. Measure the DC voltage across R4. Since the bath is requesting heat, there should be current flowing through this resistor making the voltage drop equal to 3 vdc ± 0.5 .
 - d. If the last two steps are not as stated, then most likely the CPU board is bad and it will have to be replaced.
 - e. Switch the voltmeter to an AC volts scale capable of reading 120 volts.
 - f. Measure the voltage between U1-PIN 4 and U1-PIN 6. It should be less than 1 VAC.
 - g. If it is not, then most likely the opto-coupler U1 is bad and the POWER SUPPLY board will have to be replaced. If it is, continue on.
 - h. Measure the voltage directly across the heater. It should be line voltage 110VAC to 120 VAC.

- i. If it is not, then most likely the triac Q2 is bad, and the POWER SUPPLY board will have to be replaced.
6. Check heater.
 - a. Disconnect the bath from its electrical supply.
 - b. Isolate the heater from any circuitry by disconnecting one of the heater leads.
 - c. Using an ohmmeter, check the heater resistance. Appropriate heater resistance values are listed below:

If the resistance reads 0 or infinity, then replace the heater.

	115 VOLT	230 VOLT
Model 260	15 OHMS	60 OHMS
Model 265	10 OHMS	40 OHMS
Model 270	10 OHMS	40 OHMS

- d. Check the resistance between the heater leads and ground (green wire). If the resistance reads 0 ohms, replace the heater.

***** **Problem B: Constant Heat** *****

1. Bath is heating when it is not requesting heat.
 - a. Select a setpoint temperature at least 10° below the actual temperature. The bath should not be requesting heat as indicated by the HEATER ON indicator lamp on the front panel. It should NOT be on continuously or even flashing.
 - b. Measure the voltage at U1-PIN 2 with respect to TP2-COM. It should be no less than 4 VDC.

- c. Measure the DC voltage across R4. The voltage should be approximately zero volts.
- d. If the last two steps are NOT as stated, then most likely the CPU board is bad and it will have to be replaced, otherwise continue on.
- e. Switch the voltmeter to an AC volts scale capable of reading 120 volts.
- f. Measure the voltage between U1-PIN 4 and U1-PIN 6. It should be line voltage 110 VAC to 120 VAC.
- g. If it is not, most likely the opto-coupler U1 is bad and the POWER SUPPLY board will have to be replaced. If it is, continue on.
- h. Measure the voltage directly across the heater. It should be approximately 0 VAC.
- i. If it is not, then most likely the triac Q2 is bad, and the POWER SUPPLY board will have to be replaced.

2. Check Temperature Probe. Problem A, Step 4.

*** **Problem C: Unstable Temperature Control or Display** ***

1. Use gable cover provided to improve temperature control.
2. If control is stable but not at desired temperature, then check temperature calibration. Re-adjust if necessary.
3. Check Temperature Probe, Problem A, Step 4.

**** Problem D: Display Reads NNNN or UUUU ****

The control boards have the ability to detect an OPEN or SHORTED temperature sensor. These two conditions are shown on the display as follows:

<u>DISPLAY</u>	<u>CONDITION</u>
UUUU	SHORT CIRCUIT
NNNN	OPEN CIRCUIT

Before coming to the conclusion that the probe is bad when one of these displays appear, check the connection of the probe to the circuit board for polarity and alignment.

*** Problem E: Display Reads all Eights (8888) ***

1. The most common problem cause of this problem is the failure of the driver integrated circuit on the DISPLAY/KEYBOARD board (#3176679). Replace it.

***** Problem F: Pump Does Not Circulate Water *****

1. Verify that the water level is above both the pump inlet and outlet.
2. If the pump is not running, verify that power is applied to the bath. If power is present, check wiring of pump to power switch. If wiring is correct, replace the pump.
3. If pump is running but water is not circulating, open the control panel and gently squeeze the long piece of rubber tubing a few times to help remove any air trapped within the pumping system.
4. If water is still not circulating, check for foreign matter in the pump inlet fitting.

PARTS REPLACEMENT

Before removing any parts for replacement, verify part in question by following the instructions listed in the troubleshooting guide.

Refer to Parts Replacement Table below for appropriate replacement procedures. Failure to follow parts replacement procedures may cause damage to the bath.

WARNING

BEFORE REPLACING ANY PART, TURN POWER SWITCH TO OFF, THEN UNPLUG THE UNIT FROM ITS POWER SOURCE. SERVICE SHOULD BE PERFORMED BY A QUALIFIED TECHNICIAN.

CAUTION

WHEN REPLACING TEMPERATURE SENSOR, HEATER OR CIRCULATING PUMP, MAKE SURE ALL SEALS ARE WATERTIGHT BEFORE APPLYING ELECTRICAL POWER.

PARTS REPLACEMENT	
NOTE: For easiest access to the control housing, place bath on the edge of a bench or table. Remove four screws on upper side of control housing and carefully lower hinged control panel to the vertical position.	
Circulating Pump Replacement	<ol style="list-style-type: none">1. Drain water from bath.2. Remove the two pump leads attached to power switch and green wire from pump to ground lug.3. Loosen the two visible hose clamps with a screwdriver and remove tubing from pump.4. Disconnect pump from pump bracket by removing two screws and remove pump from bath.5. Reverse above procedure to install new circulating pump.
Replace Temperature Sensor	<ol style="list-style-type: none">1. Drain water from bath.2. Remove diffuser pan from bath.3. Remove connector J302 from CPU board.4. Remove temperature sensor fitting from inside bath pan.<ul style="list-style-type: none">- Hold small nut on fitting firmly with a wrench.- Remove large fitting piece with another wrench.5. Remove rubber grommet from temperature sensor.6. Carefully slide sensor through the control housing and remove.7. Reverse above procedure to install new temperature sensor.8. Recalibrate new temperature sensor. See Temperature Calibration on page 15.

PARTS REPLACEMENT (CONT'D)	
Replace Heater	<ol style="list-style-type: none"> 1. From inside control housing, remove 2 nuts and disconnect wires from heater leads. 2. Remove two hex nuts from heater ends using a 13/16" socket and remove heater. 3. Reverse above procedure to install new heater.
Replace "HI LIMIT" Thermostat	<ol style="list-style-type: none"> 1. Remove four leads from thermostat terminals. 2. Remove two screws from thermostat bracket and remove bracket. 3. Carefully tilt bath on its side. 4. Remove eight screws from bottom plate and remove plate. 5. Remove insulation from bottom of pan and between pan and control panel wall. 6. Loosen two nuts from thermostat bulb bracket. 7. Carefully remove thermostat from control housing by sliding bulb through hole in control panel wall. 8. Reverse above procedure to install new Hi Limit thermostat. Use extreme care when installing new thermostat. A sharp bend in the thermostat capillary will crimp off flow of hydraulic fluid.
Replace Power Supply Board	<ol style="list-style-type: none"> 1. Remove connectors J101, J102, J103, and J105 from board. Use extra care when removing J102 so as not to bend connector pins. 2. Remove five nuts from board and then remove it. 3. Reverse above procedure to install new power supply board.
Replace CPU or Display/Keyboard Board See note on the following page **	<ol style="list-style-type: none"> 1. Remove connectors J302 and J303 from CPU board. 2. Remove connector J102 from power supply board. Use extreme care when removing J102 so as not to bend connector pins. 3. Remove four nuts from CPU board and remove CPU board and Keyboard/Display board. 4. Remove J201 from CPU board. Use extreme care when removing J301 so as not to bend connector pins. Remove board. 5. Reverse above procedure to install new board.

PARTS REPLACEMENT (CONT'D)																			
<p>Replace CPU or Display/Keyboard Board See note below **</p>	<p>6. The software must be configured to match the Bath model. To enter configuration mode, press the upper left pushbutton switch while applying power. The upper left switch is labeled "Set Point." This switch is used to enter configuration codes and steps through the different levels of configurations. When a triangle is shown in the chart, the pushbutton is to be depressed once for each symbol.</p> <p>Immediately following power up with the switch depressed, the readout will display a 3-digit number. This is the software revision. The readout will then begin to do a self-test of digits. Press the "Set Point" button once. The readout will display 0000. Follow the listing below to complete the configuration.</p> <table> <tr> <td>A. Enter 37 using up/down arrows.</td><td>Press "Set Point"</td></tr> <tr> <td>B. Enter 31 using up/down arrows.</td><td>Press "Set Point"</td></tr> <tr> <td>C. Enter 2.0 using up/down arrows.</td><td>Press "Set Point"</td></tr> <tr> <td>D. Enter 500 using up/down arrows.</td><td>Press "Set Point"</td></tr> <tr> <td>E. Enter 1 using up/down arrows.</td><td>Press "Set Point"</td></tr> <tr> <td>F. Display Not Used.</td><td>Press "Set Point"</td></tr> <tr> <td>G. Display Not Used.</td><td>Press "Set Point"</td></tr> <tr> <td>H. LOW TEMP CALIBRATION* LITTLE c. Factory Setting.</td><td>Press "Set Point"</td></tr> <tr> <td>I. HI TEMP CALIBRATION BIG C. Factory Setting.</td><td>Press "Set Point"</td></tr> </table> <p>Configuration is now complete</p>	A. Enter 37 using up/down arrows.	Press "Set Point"	B. Enter 31 using up/down arrows.	Press "Set Point"	C. Enter 2.0 using up/down arrows.	Press "Set Point"	D. Enter 500 using up/down arrows.	Press "Set Point"	E. Enter 1 using up/down arrows.	Press "Set Point"	F. Display Not Used.	Press "Set Point"	G. Display Not Used.	Press "Set Point"	H. LOW TEMP CALIBRATION* LITTLE c. Factory Setting.	Press "Set Point"	I. HI TEMP CALIBRATION BIG C. Factory Setting.	Press "Set Point"
A. Enter 37 using up/down arrows.	Press "Set Point"																		
B. Enter 31 using up/down arrows.	Press "Set Point"																		
C. Enter 2.0 using up/down arrows.	Press "Set Point"																		
D. Enter 500 using up/down arrows.	Press "Set Point"																		
E. Enter 1 using up/down arrows.	Press "Set Point"																		
F. Display Not Used.	Press "Set Point"																		
G. Display Not Used.	Press "Set Point"																		
H. LOW TEMP CALIBRATION* LITTLE c. Factory Setting.	Press "Set Point"																		
I. HI TEMP CALIBRATION BIG C. Factory Setting.	Press "Set Point"																		
<p>* Calibration of the low and high end of the sensor input is done at the factory. Field calibration of the sensor may be performed using the following technique: Enter the configuration routine and press the top left switch until a little "c" (Item H) is displayed. Add water that is around 10°C to the bath completely covering the sensor. Insert a calibrated thermometer. After the thermometer stabilizes, match the display setting to the thermometer reading by pressing the "UP" or "DOWN" push button switch. Press the top left switch. The readout will now display a big "C" (Item I). Remove water from the bath and add hot water of about 70 - 80°C. After the unit stabilizes, match readout display to temperature on thermometer and press the top left switch. Calibration is now complete.</p>																			

****Note:**

When replacing CPU printed circuit board reprogramming may be necessary.
Follow thoroughly instructions provided with the replacement circuit board.

TEMPERATURE CALIBRATION

The following instructions provide a step-by-step procedure for temperature calibration. ***Do not deviate from the procedure or the calibration of your bath may be inaccurate.***

CAUTION

THIS PROCEDURE SHOULD ONLY BE PERFORMED AFTER INSTALLING A NEW TEMPERATURE PROBE INTO THE BATH.

1. Connect the temperature probe connector to J302 on CPU board.
2. Turn Power switch to the OFF position.
3. Press and hold the SETPOINT key while turning Power ON.
4. The unit will now cycle through an LED display test. Press SETPOINT to end test.
5. The display will now read 0000. Press UP until display reads 0012. Press SETPOINT.

NOTE:

IF A VALUE OTHER THAN 0012 IS DISPLAYED AND SETPOINT IS PRESSED, THE UNIT WILL RETURN TO NORMAL OPERATING MODE AND DISPLAY WILL SHOW ACTUAL BATH TEMPERATURE.

6. The unit is now in calibration mode. Raw count value appears. Hit setpoint.

7. The display will now show cXX.X (XX.X is the actual bath temperature) - this is Low Temperature Calibration.
8. Place the temperature probe and a calibrated reference thermometer into an ice bath. Press UP or DOWN until the display value matches the reference thermometer. Press SETPOINT to store the new low temperature calibration value.
9. The display will now show CXX.X - this is High Temperature Calibration. (XX.X is the actual bath temperature).
10. Place the temperature probe and a calibrated reference thermometer into a steaming water bath. Press UP or DOWN until the display value matches the reference thermometer. Press SETPOINT to store the new high temperature calibration value.

NOTE


PRESSING OFFSET AT ANY TIME DURING THE CALIBRATION PROCESS WILL RETURN UNIT TO NORMAL OPERATION MODE AND DISPLAY WILL SHOW ACTUAL BATH TEMPERATURE.

REPLACEMENT PARTS LIST


	Model #260		Model #265		Model #270	
	51221035	51221037	51221039	51221040	51221036	51221038
	120V	230V	120V	230V	120V	230V
Bath Pan Kit	3167140		3167141		3167142	
Diffuser Pan Assy	3164495		3164479		33164500	
PCB Kit, PowerSupply	3166924					
PCB Assembly, CPU	3167119					
PCB Assembly, Display	3176679					
Heater Assembly Kit	3167138	3167139	3167136	3167137	3167136	3167137
Voltage Connector	3176698	3176699	3176698	3176699	3176698	3176699
Wire Harness	3176700					
Temperature Probe Kit	3167144					
'O' Ring Kit, Water Bath	3167129					
Tubing Kit	3167145					
"Hi-Limit" Thermostat Kit	3167143					
ON/OFF Switch	3175318					
Hi Limit Light, Amber	3177575	3173907	3177575	3173907	3177575	3173907
Line Cord	3178034	3176551	3178034	3176551	3178034	3176551
Pump Assembly Kit	3167333	3167342	3167333	3167342	3167341	3167343
Drain Kit (includes Nut & O-Ring)	3166191					

ACCESSORY PARTS LIST


	Model #260		Model #265		Model #270	
	3166680	3166682	3166684	3166685	3166681	3166683
Cover, Gable	3166565		3166208		3166230	
Kit, Water Level Regulator	3166223					
Sanitizer, Oakite	3166250					
Diffuser Shelf, Raised	3165861					

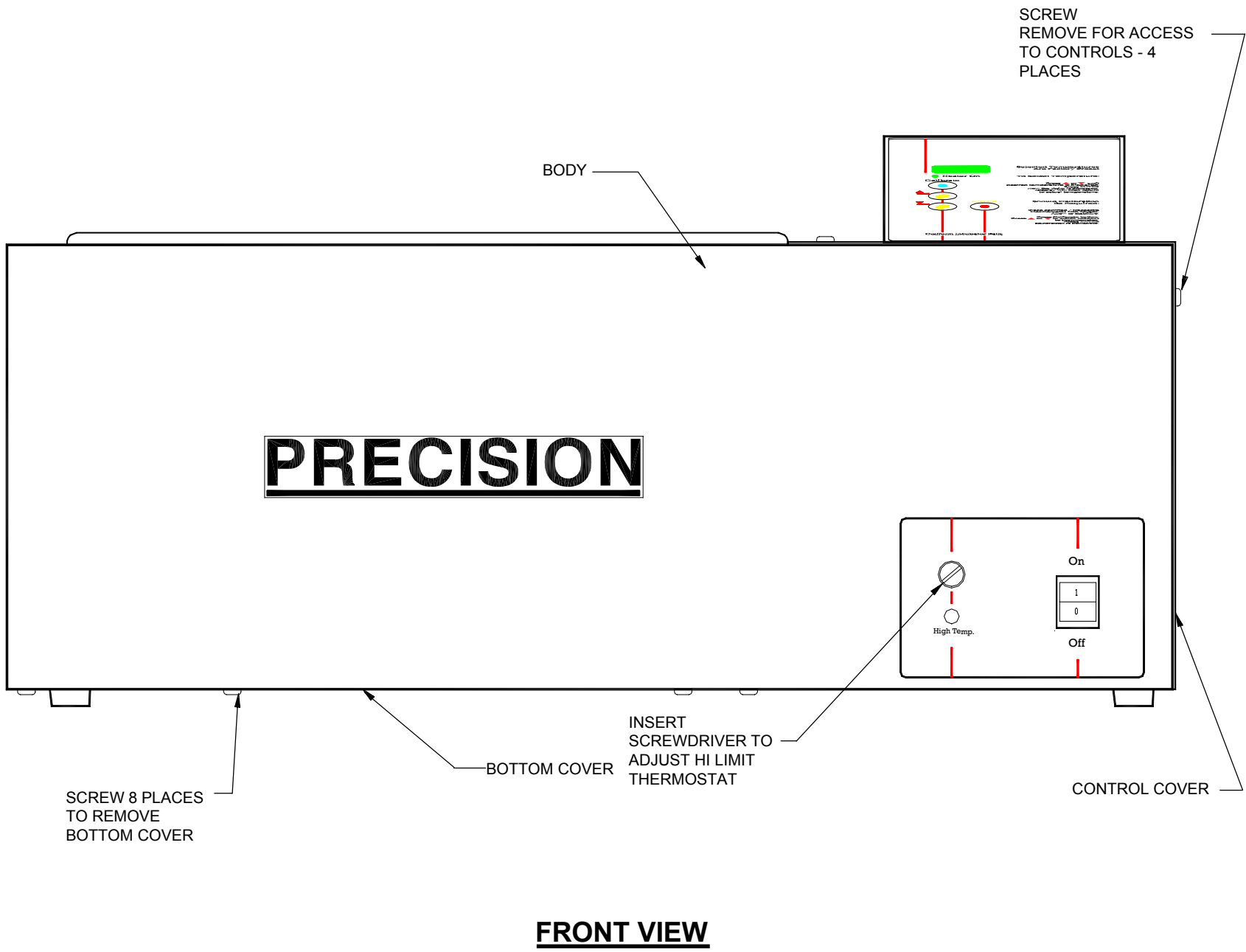


FRONT VIEW

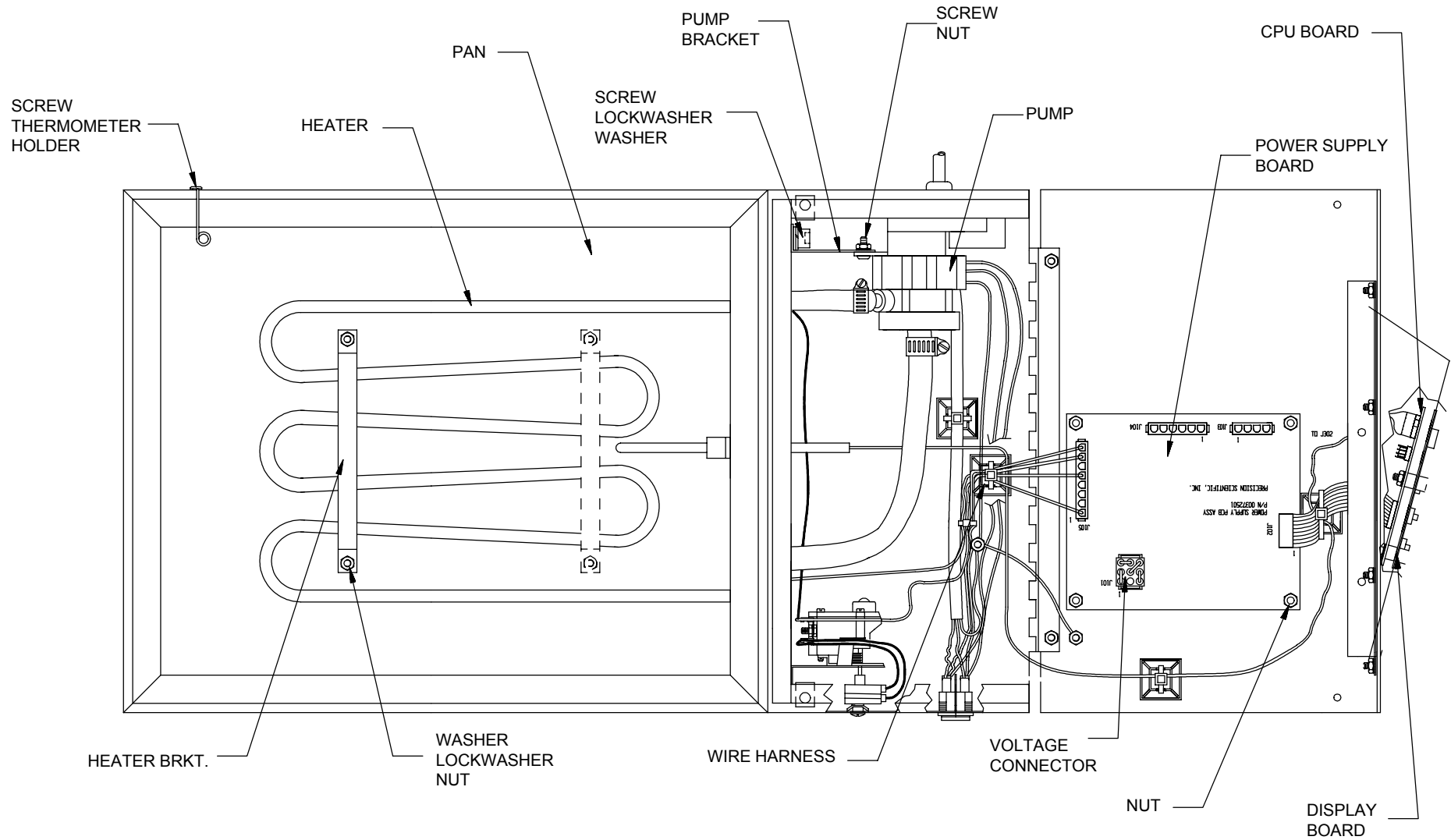


TOP VIEW
(WITH CONTROL COVER OPEN)

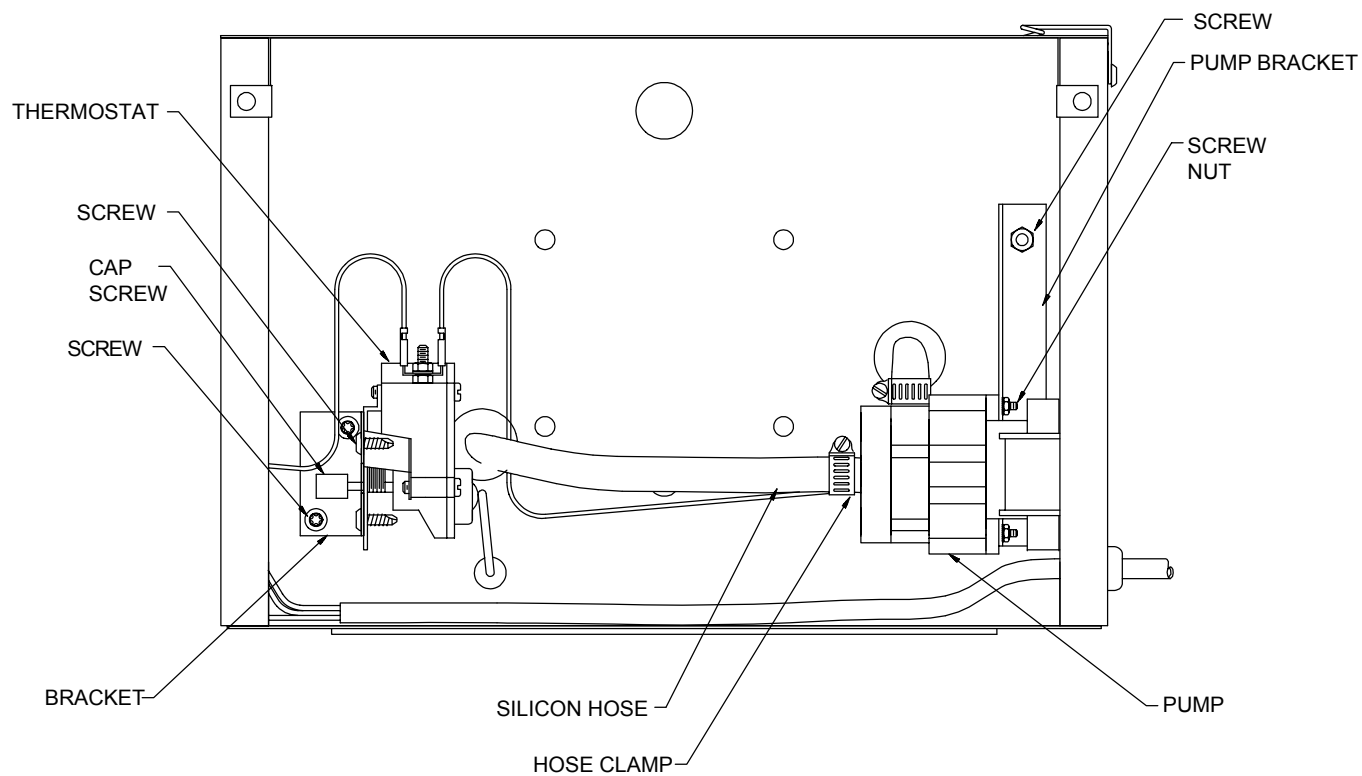


**SIDE VIEW
(WITH CONTROL COVER REMOVED AND WIRING EXCLUDED)**

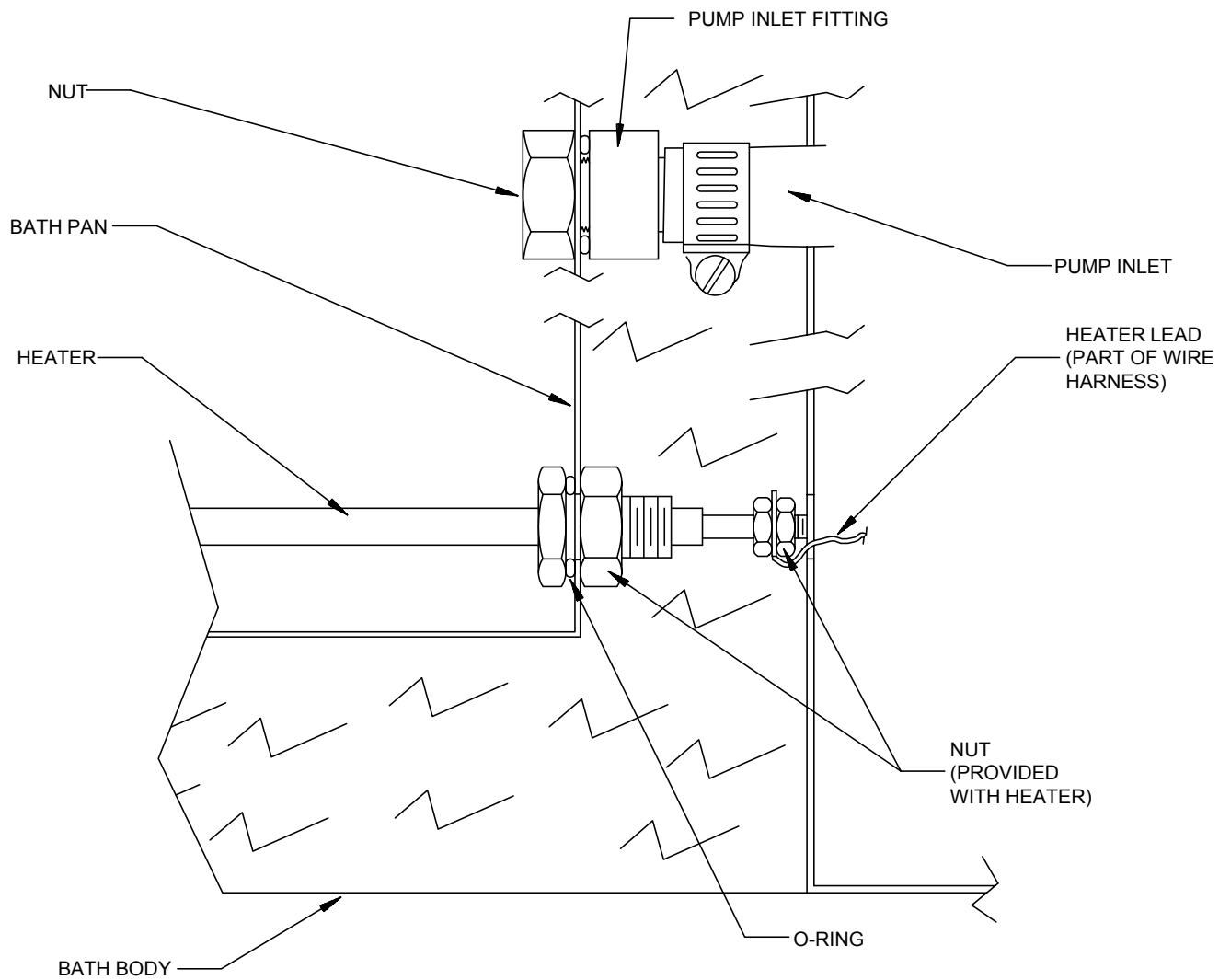


**FRONT VIEW
(CUTAWAY DETAIL OF PUMP INLET & HEATER)**

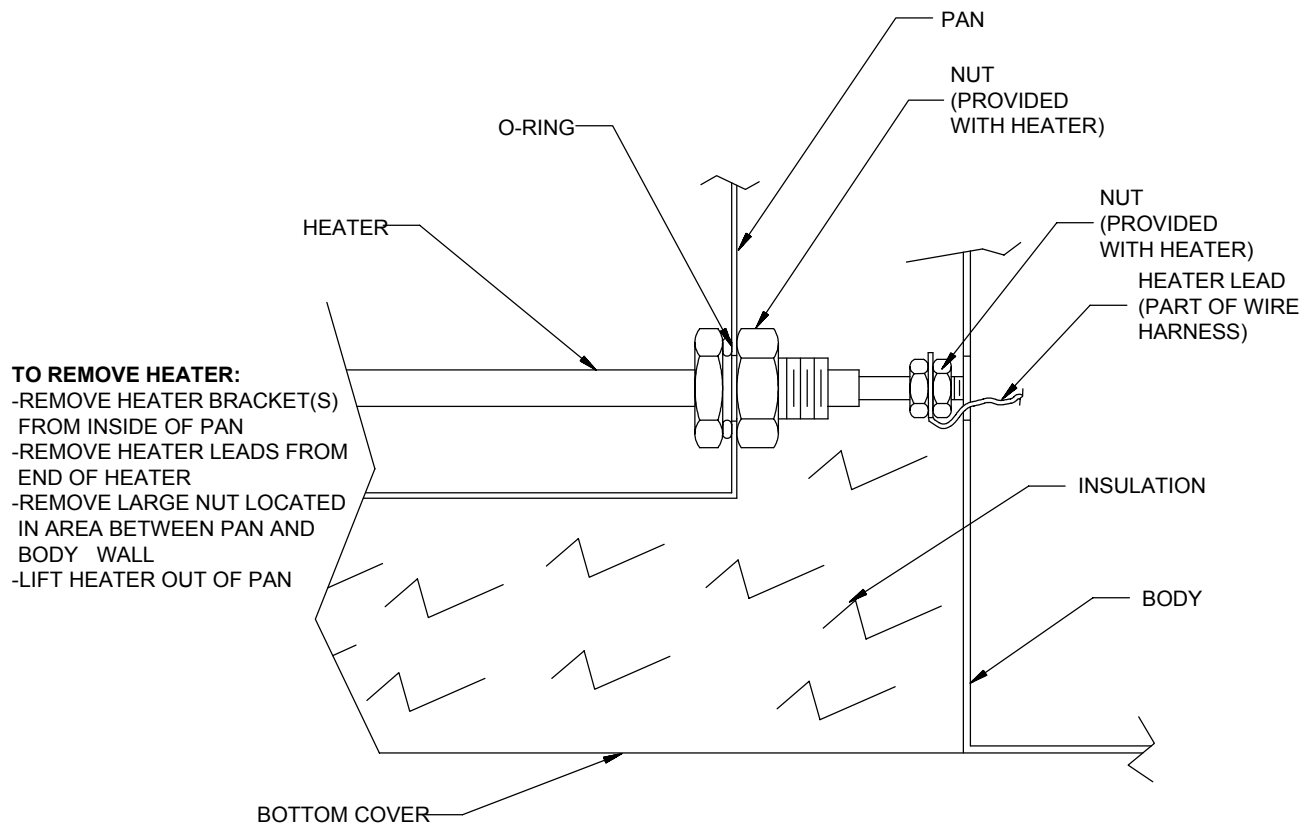


HEATER REMOVAL DETAIL

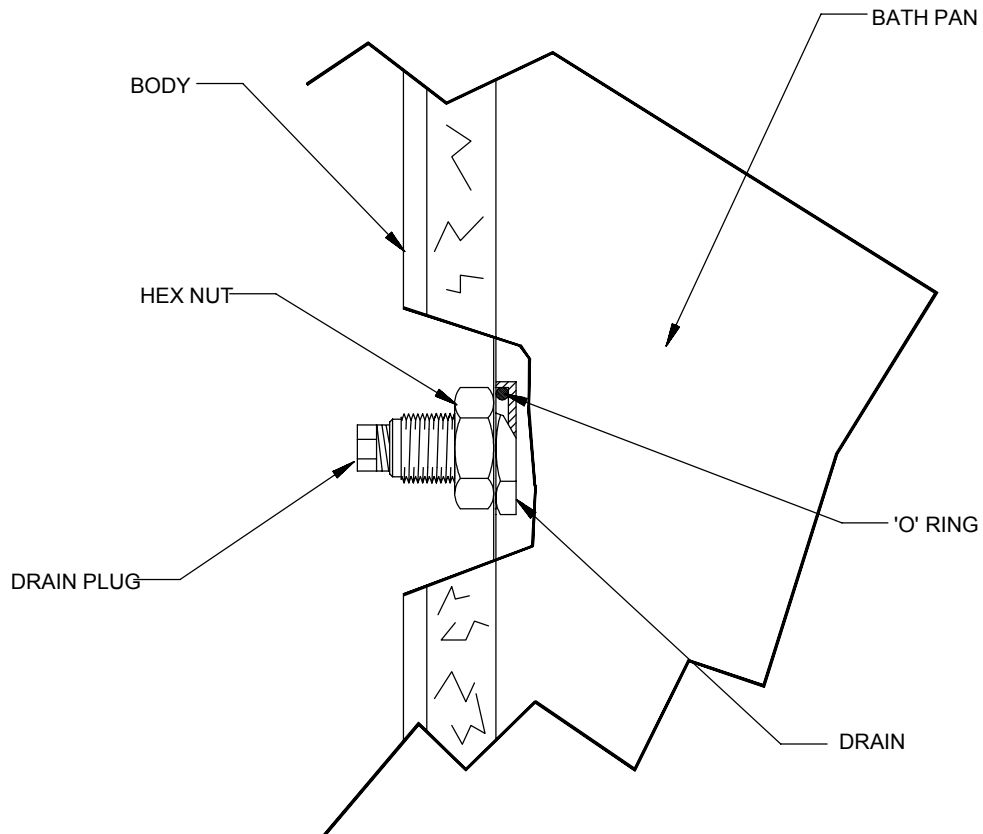


CUTAWAY VIEW OF DRAIN

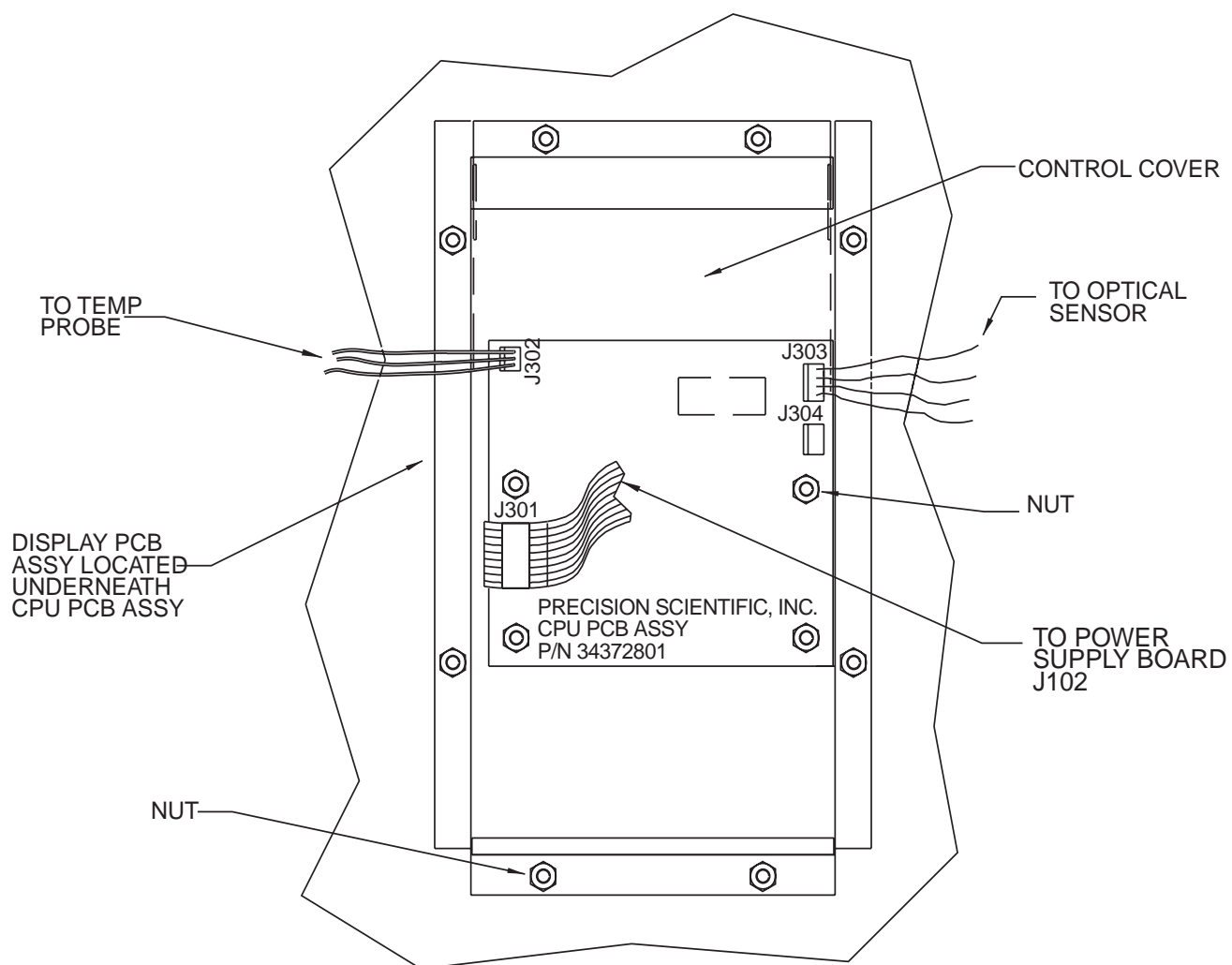

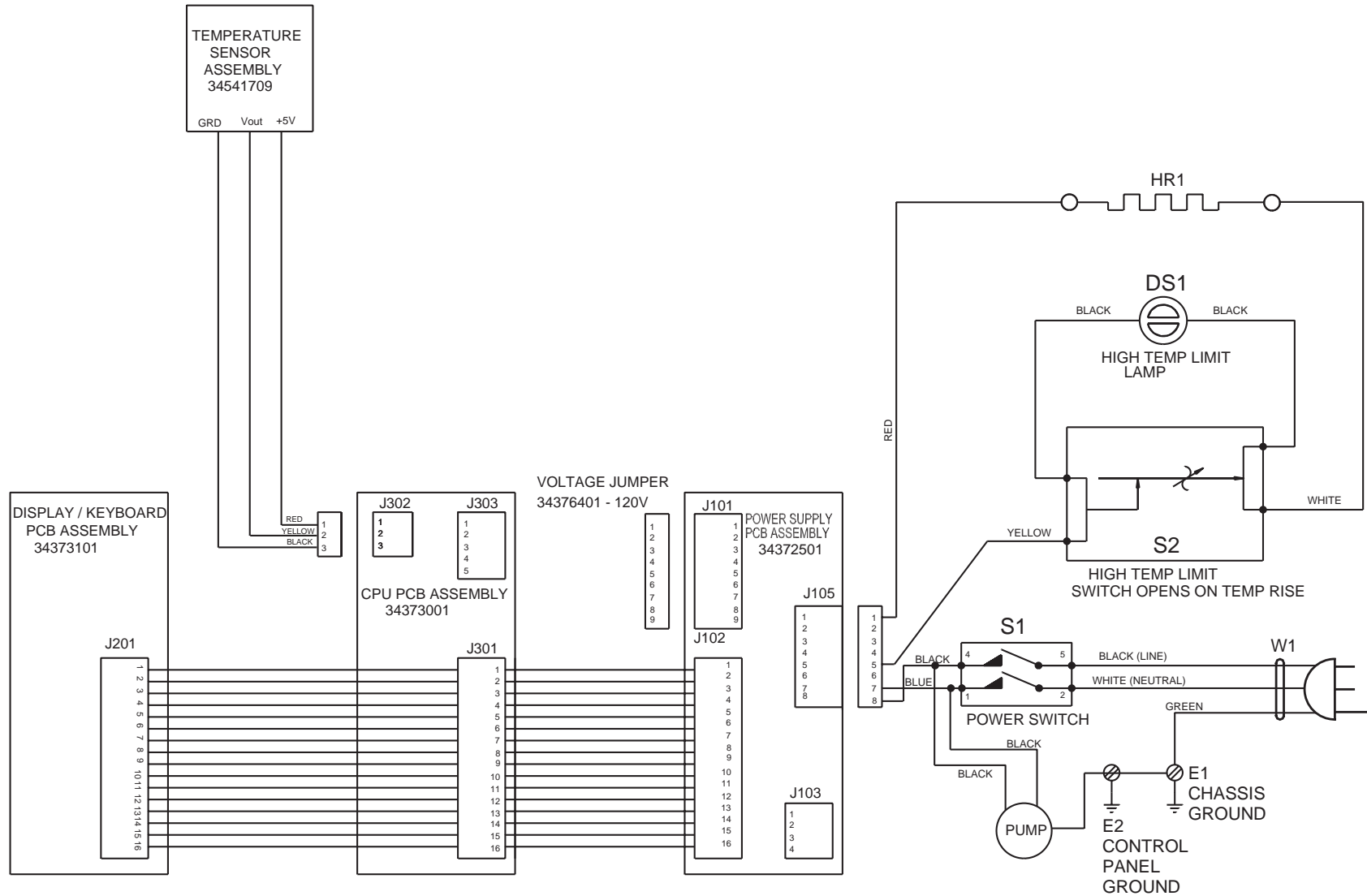

SCHEMATIC DIAGRAM
(120V)

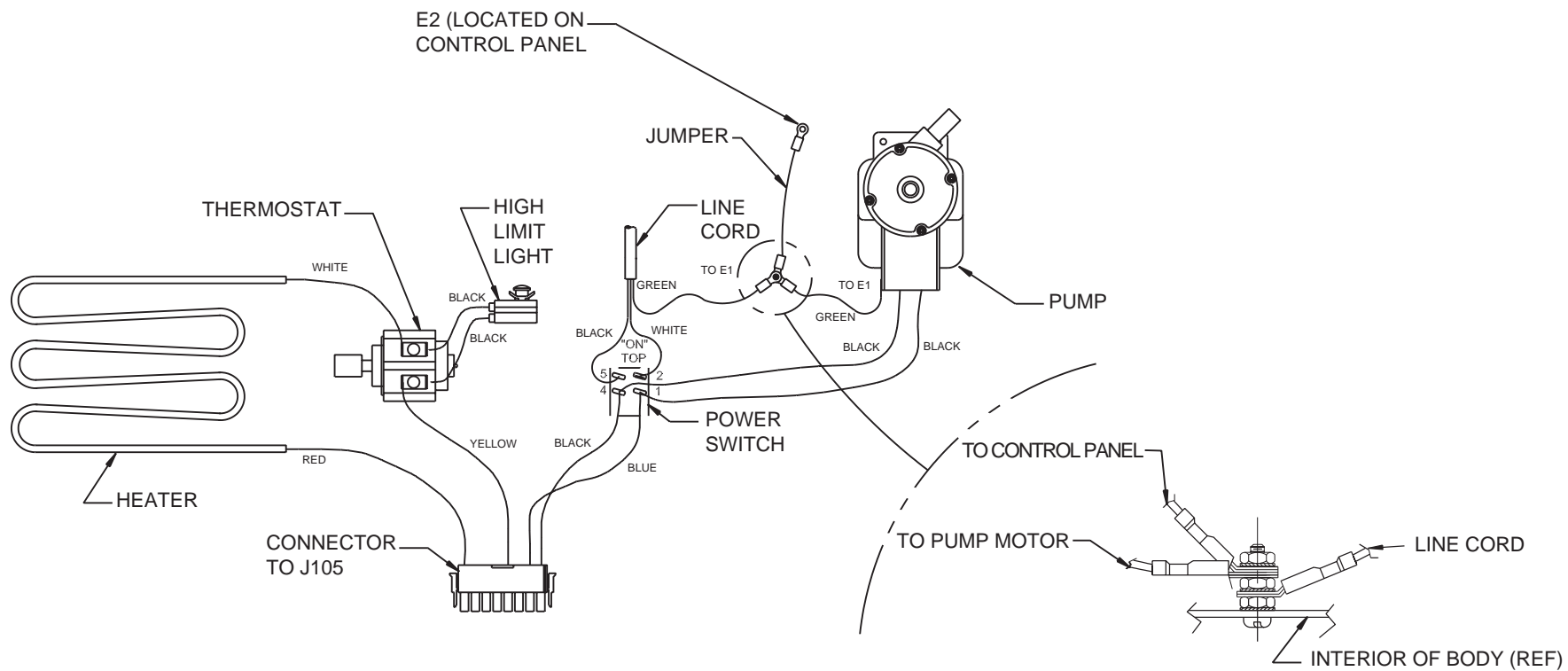

DETAIL A
GROUND CONNECTION-E2

WIRING CONNECTION DIAGRAM
(120V)s

THERMO ELECTRON CORPORATION STANDARD PRODUCT WARRANTY

The Warranty Period starts two weeks from the date your equipment is shipped from our facility. This allows for shipping time so the warranty will go into effect at approximately the same time your equipment is delivered. The warranty protection extends to any subsequent owner during the first year warranty period.

During the first year, component parts proven to be non-conforming in materials or workmanship will be repaired or replaced at Thermo's expense, labor included. Installation and calibration are not covered by this warranty agreement. The Technical Services Department must be contacted for warranty determination and direction prior to performance of any repairs. Expendable items, glass, filters and gaskets are excluded from this warranty.

Replacement or repair of components parts or equipment under this warranty shall not extend the warranty to either the equipment or to the component part beyond the original warranty period. The Technical Services Department must give prior approval for return of any components or equipment. At Thermo's option, all non-conforming parts must be returned to Thermo Electron Corporation postage paid and replacement parts are shipped FOB destination.

THIS WARRANTY IS EXCLUSIVE AND IN LIEU OF ALL OTHER WARRANTIES, WHETHER WRITTEN, ORAL OR IMPLIED. NO WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE SHALL APPLY.

Thermo shall not be liable for any indirect or consequential damages including, without limitation, damages relating to lost profits or loss of products.

Your local Thermo Sales Office is ready to help with comprehensive site preparation information before your equipment arrives. Printed instruction manuals carefully detail equipment installation, operation and preventive maintenance.

If equipment service is required, please call your Technical Services Office at 1-888-213-1790 (USA and Canada) or 1-740-373-4763. We're ready to answer your questions on equipment warranty, operation, maintenance, service and special application. Outside the USA, contact your local distributor for warranty information.



